VITEA[®] CARE

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ATTENTION

Always follow the basic precautions according to the information included below and regarding safety measures when using the device. Read the entire instructions carefully before use. Keep this user's manual for future reference.

Junior Plus children's rehabilitative stroller is available in two sizes, further known as Junior Plus 2, Junior Plus 3. For detailed information please refer to the Technical Parameters section.

1. INTRODUCTION

Thank you for choosing our product. We are deeply convinced that it will meet your expectations. Please read this document before first use. This manual was made in order to provide the necessary information on the safe use of the device. The correct application of the recommendations will contribute significantly to extended durability and visual design of the product. We would like to highlight that to davte we are constantly improving the quality of our products, which may result in minor changes that are not included in this manual.

2. PRODUCT DESCRIPTION

Fig. 1. Construction scheme of Junior Plus children's rehabilitative stroller.

1. Handle
2. Lever for changing the seat angle and backrest
3. Trunk stabilizer
4. Side stabilizer
5. Side frame
6. Rear wheel
7. Basket for small items
8. Canopy
9. Head stabilizer
10. Safety belts
11. Front guard
12. Front wheel handle
13. Footrest
14. Front wheel

EQUIPMENT OF THE STROLLER'S PACKAGE

Side frame connector

Fig 2. Stroller corpus

Fig 3. Seat

Rear wheel

Front wheel stem

Front wheel

Fig 4. Wheels

Fig 5. Front guard

Fig 6. Canopy

Fig 7. Basket for small items, umbrella (additional option), table and tray holder, lumbar support, attachable cup holder

Fig 8. Sleeping bag

Fig 9. Raincover

3. INDICATIONS

The stroller is designed for children with disabilities and / or dysfunctions movement activities as a result of damage of the CNS Central Nervous System. It can be used in children with CP Cerebral Palsy, myelomeningocele, muscle disorders, dysfunctions of the lower limbs.

The maximum weights are **40 kg** (Junior Plus 2) & **75 kg** (Junior Plus 3).

4. CONTRAINDICATIONS

Contradictions against sitting position.

5. USAGE

- After unpacking the package, check if it contains all parts in accordance with the list of elements set out above.
- Ensure that all components have no visible damage. If any doubts, please contact the manufacturer.
- Do not use the stroller until all parts are properly connected.

UNFOLDING THE STROLLER

Caregiver handle

Lever for changing angle of seat and backrest

Fig 10. Folded stroller

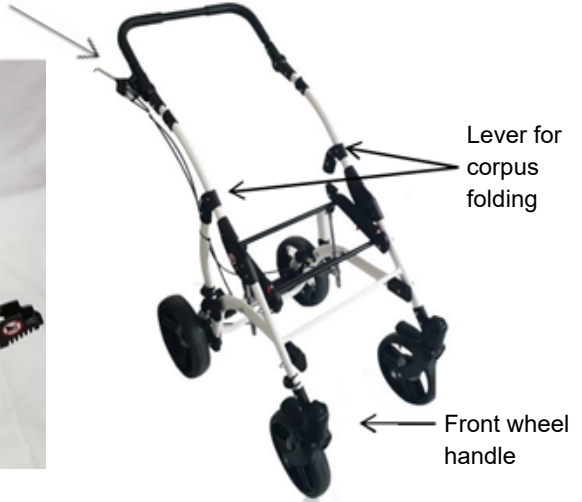
Fig 11. Unfolded stroller

- Unfolding the stroller should begin with the setting of the corpus on a flat surface, preferably on the floor. Handle to drive should be facing up (Fig. 10).
WARNING: In Junior Plus 3 start first with unfolding footrests.

- Stroller corpus should be lifted up using the handle until (Fig. 11.) you will hear a click.
- To fold a stroller corpus, at the same time pull upwards the two levers folding. At this moment lock is released and side frame connector is separated.
- Then direct down the driving handle till the stroller corpus get figure like on Fig. 10.
- Make sure that the actuators are unlocked when folding and unfolding the stroller. To do this, press the seat and backrest

WARNING! When folding the stroller, the rear wheels must be unlocked!

Fig 12. Folded stroller

Fig 13. Seat lock closeup

- Into unfolded stroller corpus, attach wheels.
- Front wheel stem slip in the wheel grip in the direction from the bottom to the top, to the point hearable click.
- To remove the front wheels, press button located at the top of wheel handle and pull it in from top to bottom (Fig. 12.).
- Rear wheels slide on the pin located in the back of the stroller corpus. Press in this same time the button located in the central part of the wheel to hearable click.
- To remove the rear wheels, press the button in the middle of the wheel and then pull the wheel out of the stroller corpus (Fig. 13.).
- Basket for small items should be mounted on the front of the stroller to a thin rod in the lower part of the side frames using latch, while at the rear of the horizontal rod using Velcro (Fig 13.).

Rod

Fig 14. Seat with rod

Seat socket

Fig 15. Seat sockets

- Seat is equipped with four pins, two on each side.
- While installing the seat, pins should be positioned in the seat socket (Fig. 15.).
- The seat should be held with both hands, one hand front edge of the seat and second hand the upper edge of the backrest.
- After placing the seat in the stroller corpus, it will lock automatically.
- To remove the seat from the corpus, the safety devices must be unlocked at the same time. The safety catch (1) needs to be moved down, and the second safety catch (2) towards the first safety catch (Fig. 16.). There will be a characteristic click. Do this on both sides. Carefully pull out the seat.

1

Fig 16. Safety catches

2

- To install the canopy push canopy handles on the stroller support frame (Fig. 17, 18.).
- Sleeping bag mount if only the child is properly fastened with the seatbelt (Fig. 19).
- To mount the sleeping bag, start from the bottom by hooking it on the footrest.
- Front guard must be attached to the seat frame with clips, make sure that the barrier is properly fastened (Fig. 18.).

Fig 17. Canopy attachment

Fig 18. Canopy folded

Fig 19. Canopy unfolded,
sleeping bag

- Raincover should be placed from top to bottom. Cover canopy and then stretch down to a footrest (Fig. 20.).
- When the raincover is no longer needed, dry it, out of the reach of children.

Fig 20. Raincover properly attached to the stroller

FOLDING THE STROLLER

- To fold the stroller, remove the canopy and front guard.
- Open seat lock, remove seat.
- Unlock the rear wheels (you can also remove the front and rear wheels).
- Pull levers folding corpus.
- Fold corpus.

WARNING! Failure to follow these guidelines can result in damage to the wheelchair.

6. SAFETY OF USE

- a) Before using the stroller, refer to this manual and keep it for future reference.
- b) Your child's safety may be compromised case of non-compliance of rules contained in this manual.
- c) Never leave a child unattended. Always watch the child when he/she is in the stroller. Caretakers are responsible for the safety of the child.
- d) The seat is not suitable for children under 6 months.
- e) Always use the seat belts to avoid serious injuries caused by the child falling out or slipping out of the stroller.
- f) It is important to fasten seat belt properly and its adjustment to the user.
- g) The seat and belts do not replace proper supervision of caretaker.
- h) While adjusting the parameters of the stroller, the child should not be in a stroller.
- i) All repairs must be carried out by qualified personnel. And during the warranty period only by the manufacturer's service.
- j) The stroller is not suitable for running or ice-skating and roller skating.
- k) The stroller is designed for transporting only one child.
- l) It is unacceptable to use the stroller, if you found visible or palpable damage, or worn out parts and stroller parts indicates a potential risk of an accident (e.g. Consumption wheels, upholstery, cracks, etc.).
- m) Always check before use if the stroller is properly unfolded. Ensure that all locks blocked before use.
- n) To avoid injury to make sure that the child is far away when the product is being folded and unfolded.
- o) When moving the stroller through doors and other passages, make sure that on both sides there is enough space to avoid injury of hands or damage the stroller.
- p) Do not allow your child to play with this product.
- q) Do not allow your child to stand on the stroller or footrest. The stroller might fall over and hurt the child.
- r) When in motion, always hold the handle of the stroller.
- s) Be especially careful while using the wheelchair near to steps and stairs.
- t) The stroller is not intended for use on stairs and escalators.
- u) Do not to move the stroller if a child is sitting in it.
- v) Do not to move the stroller by holding the front guard. If there is a need to move your stroller, you should hold the fixed elements of the frame.
- w) In order to avoid accidental rolling off the stroller, always put on the brakes and turn on the lock of driving direction when the stroller is stopped, when patient get in, get off, changing the position.
- y) The stroller is designed for children up to 40 kg (Junior Plus 2) and 75 kg (Junior Plus 3).
- z) Stroller should always be easy to fold out. In case of problems with unfolding. You should re-read the user's manual or contact the manufacturer.

- x) Do not let the child stand on the footrest. Never leave the stroller unattended on sloping ground, even if the brakes are on.
- y) The stroller is designed for children up to 40 kg (Junior Plus 2) and 75 kg (Junior Plus 3).
- z) Stroller should always be easy to fold out. In case of problems with unfolding. You should re-read the user's manual or contact the manufacturer.
- aa) Keep extreme caution while passing an obstacle. Extreme use of the stroller could cause damage or destroy it.
- bb) Never use the stroller as a seat for the transportation of people in cars, buses, or any motor vehicles.
- cc) During transport (e.g. In the trunk), do not burden stroller, it may weaken or damage its structure.
- dd) Do not leave the stroller with the user close to open sources of fire
- ee) Be aware that the surface temperature of the upholstery and the frame parts may be increased after exposure to the sun or other heat source. Furthermore, with low temperatures there is a danger of overcooling metal parts of stroller. It may cause burns or frostbite parts of the body.
- ff) Keep caution while moving on sloping surfaces. The maximum allowable slope angle up and down is 5%.
- gg) With the stroller you can only use on surfaces where all four wheels touching the ground are to use all of four wheels.
- hh) Stroller is designed for use indoors and outdoors only on paved and even surfaces.
- ii) Use of the product i.e. The icy, wet, greasy, strongly gritty surfaces is unacceptable - it creates a risk of falling.
- jj) It is extremely important to keep all strollers' components clean in order to maintain the product parameters.
- kk) Do not add additional loads (bags, etc.) on the stroller handles.
- ll) Storing the stroller in areas with very high humidity (baths, laundry, etc.) or in the open air, also the lack of hygiene causes rapid loss of operating parameters and aesthetic of the product.
- mm) Do not leave the child in the stroller for longer, because persons with neurological affections who stay in unchanged positions may experience problems with skin nourishment and develop bedsores.
- nn) When the footrests will not be removed before folding the stroller, it can cause damage.
- oo) During use, handing of the product and during its assembly, regulating mechanisms, it might happen trapping and / or squeeze user's or caretaker's body part in holes / gaps between the elements.
- pp) Perform these steps with extreme caution. After adjustment, stabilize the position by carefully tightening the nuts / bolts.

7. EVERYDAY USE AND MAINTENANCE

TWISTING MOTION LOCK - DRIVING ONLY IN STRAIGHT

- The front wheels are equipped with twisting motion lock (Fig. 21.).
- To turn it on, turn the lock to the right on both front wheels.
- Then set the wheels in the straight-ahead, it will block the twisting motion.
- Make sure that the lock is working properly.
- To turn off the twisting motion lock, spin the lock on both wheels to the left.

Twisting motion lock - turn on

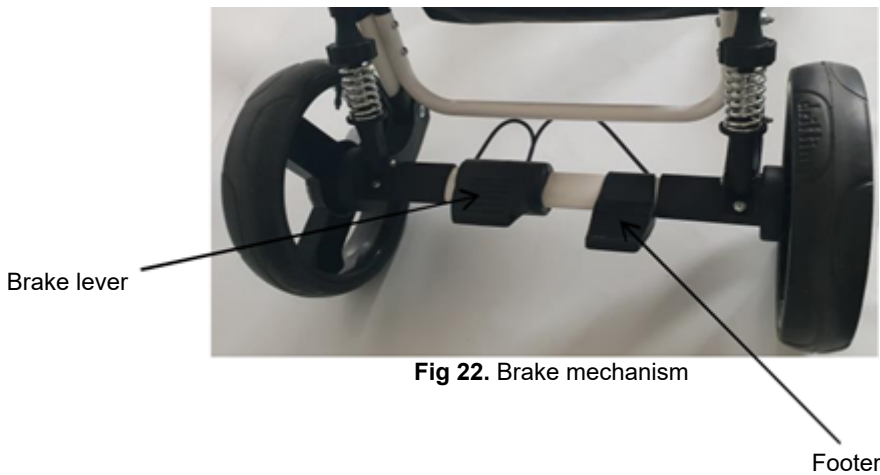
Twisting motion lock - turn off

Fig 21. Motion mechanism

SWITCHING THE BRAKE

- The rear wheels are equipped with brake system.
- To prevent unplanned movement of the stroller, switch on brakes on two rear wheels (Fig. 22.).

- To prevent movement, press the foot brake lever.
- It is necessary to make sure that the brake lever has been launched on both wheels.
- To unlock the brake, lift the brake lever.

FOOTREST LENGTH REGULATION

- The footrest has a five-stage (or seven stage) length adjustment.
- To change the length of the footrest, remove two footrest locks on both sides of the footrest (Fig. 23.)
- Select required length on the footrest, by using the holes in footrest frame.
- Next insert the footrest locks into the holes.
- It is important that each side footrest locks were placed in holes at the same height.

Fig 23. Footrest regulation scheme

Footrest locks closeup

- Select required length on the footrest, by using the holes in footrest frame.
- Next insert the footrest locks into the holes.
- It is important that each side footrest locks were placed in holes at the same height.

FOOTREST POSITION ADJUSTMENT

- By pressing by both sides of the footrest adjustment button, you can change its position up or down (Fig. 24.).
- It is important that the regulation buttons are pressed by thumbs, while the rest of the hand, grasps the footrest frame.

Fig 23. Footrest regulation scheme

SEAT DEPTH ADJUSTMENT

- Six knobs are located under the seat.
- Their use has been shown on figure 25.
- To adjust the seat depth, loosen the two knobs the seat depth adjustable.
- Grab the seat edge by hand, push or pull to set the desired depth.
- Then tighten the two knobs to adjust the seat depth to prevent automatic changes the depth of the seat.

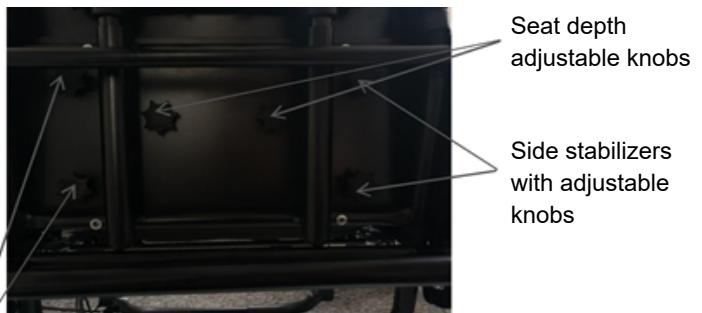
Fig 24. Seat depth adjustment

Side stabilizers with adjustable knobs

FOOTREST LENGTH REGULATION

- After unzipping in the rear of the backrest is visible eight knobs. Their use has been shown in figure 26.
- To change the position of the side stabilizers, should loosen the right pair and left pair of side stabilizers adjustable knob.
- Grasp side stabilizer on the left and then on the right, by using hand, set the desired position of the stabilizers.
- Tighten both knobs which were loose before.

FOOTREST LENGTH REGULATION

- After unzipping in the rear of the backrest is visible eight knobs. Their use has been shown on figure 26.
- To change the position of the trunk stabilizers, should loosen the right pair and left pair of trunk stabilizers adjustable knob.
- Grasp trunk stabilizer on left and then on right by using hand, set the desired position of the stabilizers.
- Tighten both knobs which were loose before

Fig 26. Motion mechanism

BACKREST HEIGHT ADJUSTMENT

- To change the backrest height, should loosen pair of backrest high adjustable knob (Fig. 26.).
- Grasp upper edge of the backrest by using hand, set the desired position.
- Tighten a pair of knobs which were loose before.

HEIGHT ARMREST ADJUSTMENT

- To change the height of the armrest, unscrew both knobs located under the armrest, then select the height and tighten the knobs (Fig 27.).

Armrest height adjustable knob

Fig 27. Armrest adjustment

ADJUSTING THE BACKREST ANGLE

- To change the angle of the backrest, push by thumb on the upper edge of the backrest. Grab lever adjusting the angle of the backrest from the bottom. Pull the lever and set convenient location in one of four positions.

SAFETY BELT SYSTEM

Fig 28. 5-point safety belts

- The stroller is equipped with a stabilizing belt which length can be adjusted.
- Settings of all stabilizing elements should be discussed with your doctor or physiotherapist.
- To ensure safety, after placing the child in the stroller, caretaker should immediately buckle them into stabilizing belts.
- Security system include three-point seat belts with central buckle and lap belt with buckle.
- All elements of the buckle should be placed in the central fastener. To do this, press the button on the buckle and place the elements (Fig 28, 29.).
- While placing elements of buckle into central buckle there must be a click. Only if click is heard, 5-point belts are properly fastened.
- To properly protect your child, adjust the length of the belt using the adjustment belt elements.
- To unfasten the seat belt, press the button on the central fastener.

Fig 28. 5-point safety belts

ADJUSTING ANGLE OF THE SEAT AND THE BACKREST

- The stroller is equipped with a lever to change the angle of the seat and backrest (Fig 29A.).
- To change the angle of the seat and backrest, push by left hand the lever change the angle of the seat and backrest into the stroller driving handle.
- By right hand grab the upper part of the backrest and move up or down to set the desired setting.

Height adjustment button handle

Lever for changing angle of the seat and backrest in **Junior 2 & 3**

Fig 29A. Angle adjustment

- **Junior Plus 3** is equipped with the additional mechanism that allows to adjust angle of the backrest (Fig 29B.).

Fig 29B. Angle adjustment

DIRECTION OF TRAVEL ADJUSTMENTS

- Stroller offers the ability to change settings for the seat relative to the direction of travel.
- To make changes, open the seat lock on both sides of the stroller.
- Remove the seat from the body corpus. Turn the seat by 180 ° relative to the corpus.
- Place the seat on the corpus in such a way that the „rods” are placed in the seat sockets (Fig 30). To make changes, open the seat lock on both sides of the stroller.
- Remove the seat from the body corpus.
- Turn the seat by 180 ° relative to the corpus.
- Place the seat on the corpus in such a way that the „rods” are placed in the seat sockets (Fig 31.).

Fig 30. Rods installation

Fig 31. Rods installation

TRAY INSTALLATION

- To assemble the tray, remove the front guard and then mount the table holder in the same place.
- Insert the tray into the holes in the handle.
- When dismantling the tray, press the side part of the tray in the hole from the side. Do this on both sides of the handle.
- Remove the handle.
- The tray distance can be adjusted using two screws located on the inside of the handle. Unscrew both knobs, pull out the screws, then move the part with the hole to the tray to the desired depth, insert the screws and tighten the knobs (Fig 32, 33.).
- Do this on both sides. The depth must be the same.

Fig 32. Holes for the table

Fig 32. Properly attached try

WARNING! When changing the angle of the seat and backrest, be careful of items on the tray. Too big angle change can cause objects to fall. The option with a table is not possible when you change the seat settings relative to the direction of travel.

MOVING ON SLOPES

- While moving on slopes, use caution and reduce your speed.
- Do not move on too high slopes. The maximum allowable slope angle up and down is 5%.

ENTRANCE TO THE THRESHOLD OF CURBS

- To enter stroller to the threshold or curb, it must be set so that the front wheels were set on straight ahead threshold or curb, and near to it.
- Then place a foot on the lock back of stroller, press lightly and push the handle down .
- The front wheels of stroller will lift slightly upwards, what allows to overcome the obstacles. Once the front wheels are on the obstacle, raise the rear wheels by tilting the stroller forward.

PUMPING TYRES

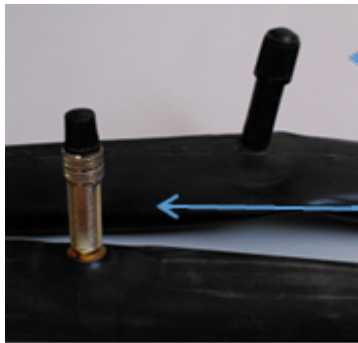
Fig 34. Valves with cover

Schrader valve

Presta valve

Fig 35. Valves w/o cover

- Before pumping the tires, identify the type of valve which is located on the tire. This may be one of the two valves shown in Illustration 34. and 35.
- Presta valve also has an additional small top, which must be loosened before pumping.
- After identifying the valve, prepare the pump designed to valve type
- Most of the new pumps are equipped to service both Presta valve and Schrader. Many of them let to reverse the tip. Smaller hole is Presta valve, Schrader valve greater. If the pump does not support Presta, you will need to buy for yourself the „adapter Presta”, which is available at any bike store.
- At this point, if the wheel is provided with a valve Presta, it will be necessary to unscrew the top located on the vent. It is a small top, which can be unscrewed with fingers. When doing this, the top should move upwards. Then valve is considered open.
- Then remove the plastic top from the valve and set it aside in a safe place. It is important to not lose it, because the valves will get dirt, leading to difficulties in pumping and rapid air loss.
- Place the nozzle of the pump to the valve and check if it stick properly. If it is needed, push the pump to the valve to make sure that the air does not escape to the sides.
- If the pump has a lever with clamp at the nozzle, tighten it. If the nozzle of the pump and valve have a threaded tip, screw them together. As a result, air will be streamed directly to the inner tube and the nozzle cannot detach from the valve.
- At the beginning of pumping the pressure in tyres should be checked at some points. Check tire pressure with your thumbs at the end of the tread. If it is possible, make only a millimetre indentation. It means that the tire has enough air. If necessary, continue pumping.
- After pumping, remove or unscrew the nozzle from the valve outlet while compressing it is ends by thumb. It is recommended to do this quick enough in order to not lose too much air. For Presta valve, tighten previously loosened additional small top. Screw the plastic top on the valve.

MAINTENANCE AND CONSERVATION

- Clean the plastic parts and the frame with damp cloth without the addition of detergents or other similar ones.
- Systematically remove the dirt from the movable elements of the stroller (screws, resistance mechanisms). Regularly remove dirt from the upholstery with a damp cloth.
- Keep clean all strollers' elements for to save functional parameters of the product.
- Storing in rooms of very high humidity (baths, laundries etc.) or in the open air, in the combination with the lack of hygiene will cause the fast loss of functional and aesthetic parameters. Do not expose the wheelchair on to extreme weather conditions (rain, snow very strong sun), because certain parts can change color.
- The stroller has coverage (coverage is part of the upholstery having the closest contact with the child, does not have stiffeners), which can be washed by hand at a 30°.

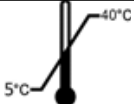
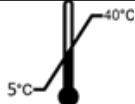
You must not use:

- solvents,
- toilet detergents,
- sharp brushes and hard objects,
- detergents that contain chlorine,
- corrosive detergents,
- clean the stroller with pressure or steam aggregates.

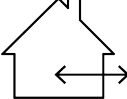
Frequency	Object of control	Comments
Before first installation	Checking the technical condition and all the elements within the packaging.	Check if any elements have any visible mechanical damage.
Before each use	Control of wheels, brakes, the lock of driving direction	The brakes can wear out, check the condition of the brakes and the lock of driving direction before each use You should also check the condition of wheels and clean them from dust and sand.
Every two weeks	Screw and nut control regarding their correct fixing.	Make sure if all the screws and top are fixed tightly.
Every two months	Cleanliness and general condition.	In case of mechanical damage contact the producer immediately.
Once a month	wheel axles	Wheel axles must be conserved, oil between the axle and the hub a few drops of lubricant.

8. TECHNICAL PARAMETERS

Parameter	Junior Plus 2	Junior Plus 3
Permitted user weight	30 kg	75 kg
Maximum load of basket for small items	1 kg	1 kg
Maximum width of the stroller	60 cm	61 cm
Maximum length of the stroller	100 cm	105 cm
Maximum height of the stroller	116 cm	120 cm
Dimensions of the stroller folded with wheels	100x60x116 cm	105x61x120 cm
Seat width	26-37 cm	32-43 cm
Seat depth	29-35 cm	34-40 cm
Backrest height	50-55 cm	60-70 cm
Front wheels diameter	22 cm	22 cm / other
Rear wheels diameter	27 cm	27 cm / other
Weight of the stroller	19,5 kg	21 kg
Maximum patient height	130 cm	165 cm
Footrest	footrest with 5-step height adjustment (every 25 mm), with foot straps and shin straps from 24 cm to 34 cm	footrest with 7-step height adjustment (every 25 mm), with foot straps and shin straps from 26 cm to 40 cm
Armrest height	18-23 cm	22-27 cm
Type of wheels	PP/PU	PP/PU
Material	Aluminum, steel, plastic, rubber	Aluminum, steel, plastic, rubber
Maximum safe slope		
Description of the environment in which the product will be used		
Temperature of storage and use		
The humidity for the storage and use	30% - 70%	30% - 70%
Terms of washing upholstery		
Serial number	Given on the product	Given on the product

9. SYMBOLS

	Confirmation of product compliance with EU standards
	Production date
	Manufacturer
	See the instructions for use
	WARNING
	Catalogue number
	Serial number
	Protect from moisture
	Protect from sunlight
	Permissible user weight
	indoor and outdoor use

	Possibility of fingers jamming
	Temperature of storage and use
	Maximum safe slope
	The product is intended to transport by car
	The product is not designed for use as a seat to transport in a motor vehicle
	It is recommended to do a short and careful hand wash at 30 ° C, do not use bleach, do not iron, do not dry clean, do not use stain removers containing solvent, not use tumble dryer
	Unique Device Identifier
	Medical device

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mdh Co., Ltd.